



**Internal Storage:**

<b>1200</b> NIBP Measurement Review	<b>120 h</b> Trend Review	<b>60</b> Alarm Review	<b>120 s</b> Frozen Waveform
--	------------------------------	---------------------------	---------------------------------



www.mennenmedical.com

Israel: Mennen Medical Ltd. 4 Ha-Yarden Street, Yavne, P.O. Box 102, Rehovot 7610002, Israel. Phone: +972-8-9323333. Fax: +972-8-9328510  
 United States: Mennen Medical Corp. 290 Andrews Road, Feasterville-Trevose, PA 19053-3480, USA. Phone: +215 259-1020. Fax: +215 357-2010  
 United Kingdom: Charter Kontron Ltd. Unit 18 Avant Business Center 21 Denbigh Road Milton Keynes MK1 1DT, UK. Phone: +44 1908 646070. Fax: +44 1908 646030



We Follow Every Beat of Your Heart

**MennMove 50**

Portable Patient Monitor with easy remote access




**8.4" Color TFT-LCD Display**

**Standard Parameters**  
3/5-lead ECG, HR, RESP, SpO2, NIBP, PR, 2-TEMP

**Optional Parameters**  
Nellcor OxiMax SpO2, 2-IBP, Respironics CO2, Touchscreen, Thermal Recorder, W-LAN Accessories Kit

**MennMove 70**

Patient Monitor

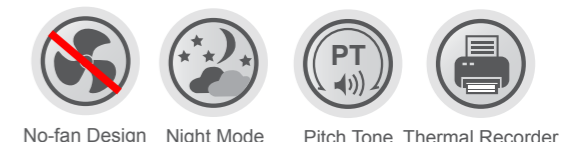

**12.1" Color TFT-LCD Display**

**Standard Parameters**  
3/5-lead ECG, HR, RESP, SpO2, NIBP, PR, 2-TEMP

**Optional Parameters**  
2-IBP, C.O., CO2 (Sidestream), Respironics CO2 (Mainstream and Sidestream), PHASEIN Anesthesia Gas/O2, Nellcor OxiMax SpO2




# MennMove 50



## Reliable Algorithms

- iSEAP™ ECG algorithm optimized for arrhythmia detection, pacemaker detection, and HR measurement
- iMAT™ SpO2 algorithm with outstanding motion resistance and low perfusion resistance performance
- iCUFFS™ NIBP algorithm optimized for cardiac patients, hypertensive patients, and neonatal patients

## Multiple Display Modes

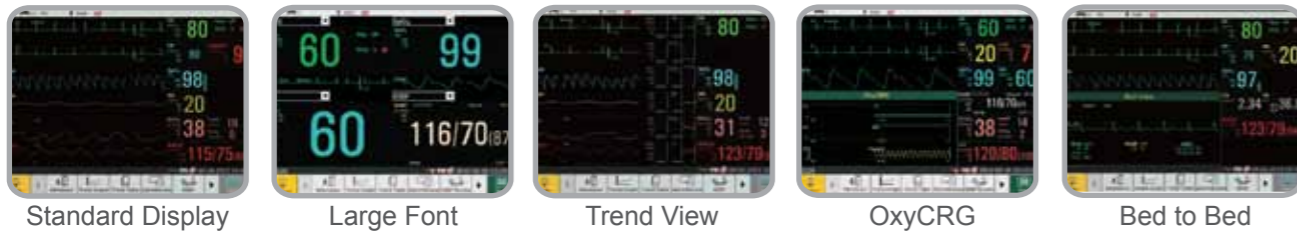

Standard Display    Large Font    Trend View    OxyCRG    Bed to Bed

## 3G Ambulatory Care




Ensemble CMS

Web Viewer

## Respironics CO<sub>2</sub> (Mainstream/Sidestream) for Intubated/Non-intubated Patients



- Plug & Play
- Dehumidification tube instead of water trap design
- No need to calibrate on a regular basis
- Low sampling rate of 50 ml/min

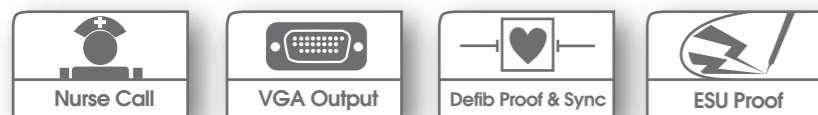

# MennMove 70



## Reliable Algorithms

- iSEAP™ ECG algorithm optimized for arrhythmia detection, pacemaker detection, and HR measurement
- iMAT™ SpO2 algorithm with outstanding motion resistance and low perfusion resistance performance
- iCUFFS™ NIBP algorithm optimized for cardiac patients, hypertensive patients, and neonatal patients

## Multiple Display Modes



Standard Display    Large Font    Trend View    OxyCRG    Bed to Bed

## End tidal Carbon Dioxide for Intubated/Non-Intubated Patients



### CO<sub>2</sub> (Sidestream)

- Superior water trap design for accurate monitoring
- iCARB™ algorithm with intelligent CO<sub>2</sub> pseudo wave identification technology
- Multiple sampling accessories as options for adult, child and neonate patients

### Respironics CO<sub>2</sub> (Mainstream/Sidestream)

- Plug & play module design
- Dehumidification tube instead of water trap
- Low sampling rate of 50ml/min suitable for all types of patients

## PHASEIN Anesthetic Gas/O<sub>2</sub> (Mainstream/Sidestream)



- Unique Plug & Play module for mainstream anesthetic gas monitoring
- The no moisture sampling line design suitable for all clinical applications
- A stable low flow of 50mL/min for all types of patients

