

Monitor Configuration and Options:

- A multi-parameter vital signs monitor
- ECG, Respiration, SpO2, NIBP, Temperature
- Optional Invasive blood pressure (2 or 4) channels and thermo-dilution cardiac output
- Optional Built-in recorder
- Touch screen
- Optional Microstream EtCO2
- Line operated with 3 hours (optional 6 hours) battery backup
- Data transfer to central station and printer via wired or wireless (optional) Local Area Network
- Long term memory as shown in the following table:


Storage Type	8GB	16GB	32GB	64GB
Charts	80 Days	80 Days	80 Days	80 Days
Trends	80 Days	80 Days	80 Days	80 Days
Full disclosure	~5 Days	~8 Days	~8 Days	~8 Days
Overview	~5 Days	~8 Days	~8 Days	~8 Days
Event strips	75	75	75	75

Features & Benefits:

- Long-term data storage:
 - ◊ Graphic Trend
 - ◊ Numeric Charts
 - ◊ ECG & Vital Signs Waveform Full Disclosure
- LAN connectivity / Wireless LAN Option
- Interface to Nurse Central Station
- Full Connectivity with Hospital network
- Easily customizable for optimal clinical use
- User friendly, intuitive interface, including easy access "Quick Keys"
- **Special features** - CRG, Alarm Watch, CerebraLogik line - Dual Channel EEG & aEEG
- **Multiple Interface Options** - Ventilators, Anesthesia machine, PiCCO, Vigilance/Vigileo monitors, Radical 7, BIS
- **Customizing Tools** - Monitor Profile, "Quick keys", Permission Editor



Big numbers display

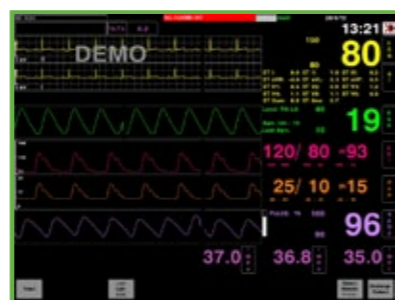

Standard display



www.mennenmedical.com

Israel: Mennen Medical Ltd.
6 Ha-Kishon Street, Yavne.
P.O. Box 102, Rehovot 7610002, Israel.
Phone: +972-8-9323333
Fax: +972-8-9328510

United States: Mennen Medical Corp.
290 Andrews Road,
Feasterville-Trevoze, PA 19053-3480, USA.
Phone: +215-259-1020
Fax: +215-357-2010

United Kingdom: Charter Kontron Ltd.
Unit 18 Avant Business Center 21
Denbigh Road Milton Keynes MK1 1DT, UK.
Phone: +44-1908-646070
Fax: +44-1908-646030




VitaLogik 4000

Compact Mobile Monitor



VitaLogik 6000

Configured Patient Monitor

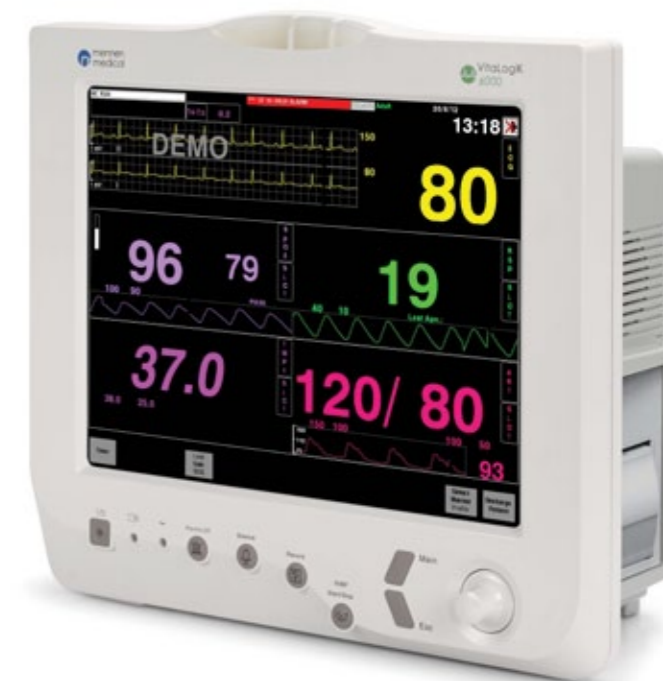

Applications:

Working successfully in all Hospital units, including:

- Emergency Room
- Trauma Unit
- Operating Room & Recovery
- Neonatal ICU
- Respiratory ICU
- Cardiac ICU

These products are CE certified
FDA Cleared


Mennen Medical Ltd. All rights reserved. Data is subject to be changed without notice.

CerebraLogik (aEEG)


The CerebraLogik is a stand alone, dual channel differential EEG amplifier, that interfaces with **VitaLogik** monitors.

CerebraLogik interface cable is connected to the UIM input of the monitor, and can be activated via the UIM menu.

The **VitaLogik** monitors display real time dual channel EEG and long term history of the Amplitude integrated EEG (aEEG). It shows aEEG on a 3 hour panel. The EEG and aEEG can be monitored and stored simultaneously with the other vital signs monitored by the **VitaLogik** monitor.

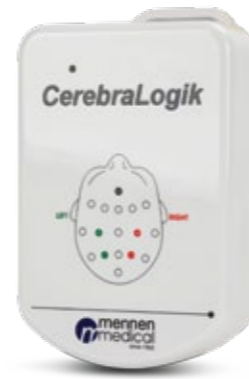

VitaLogik 4000 / 6000 Display Options & Patient Data Review

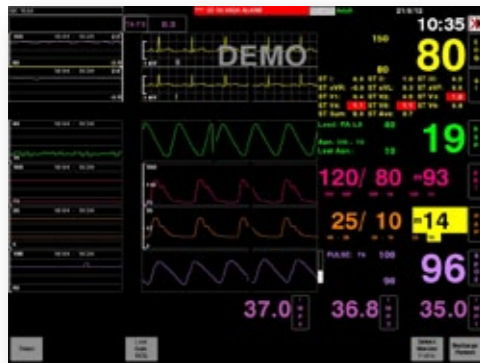

Waveform + Trend Display



Remote Bed View and Alarm Watch




Full Disclosure of ECG waveforms
10 second display



Overview of Vital Signs waveforms
synchronized to top displayed ECG



Vital signs chart with up to
one minute resolution

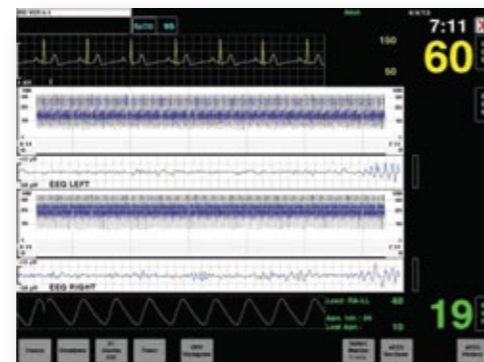

CerebraLogik aEEG Display

VitaLogik 4000



VitaLogik 6000

